



Introducing GammaVision Report Writer V3.10

Robust GammaVision Analysis. . .

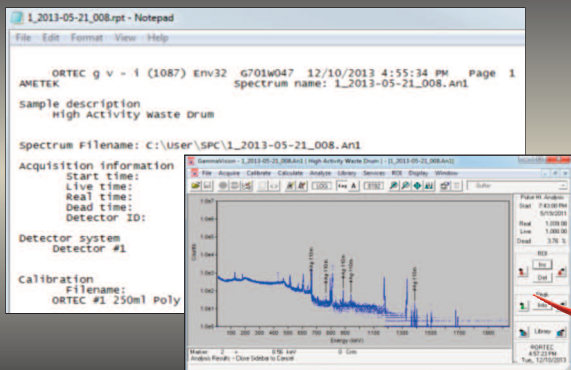

. . . in Brilliant Custom Reports

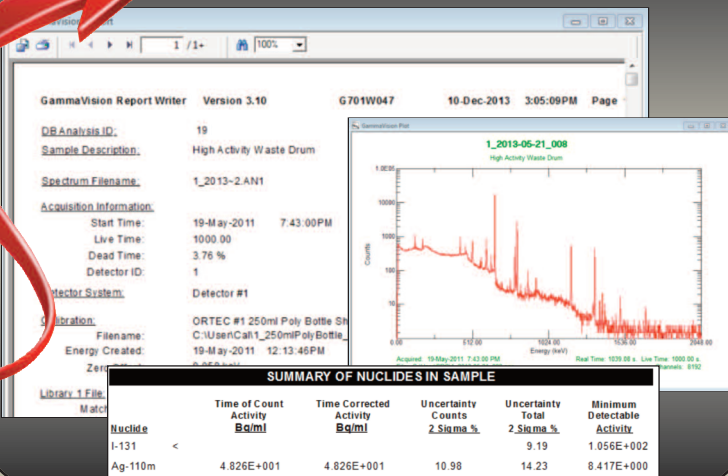

Export to Word, Excel, PDF, and More!

Why Report Writer?

Compatibility


- Windows 7 (32 & 64 bit) and Windows XP
- Microsoft Access Database
- Crystal Reports 11.5
- GammaVision 8 and earlier

Process Efficiency

- Automated Custom Reports with Spectrum Image
- Extensive Report Export Options (Word, Excel, PDF)
- Simplified LIMS Data Integration
- Data Management for Search, Archive, Purge functions

Defendable Results

- Populated directly from GammaVision Analysis
- Automated Calculations and Data Validation through Custom Report Templates
- Analysis Data supported by Spectrum Graphics

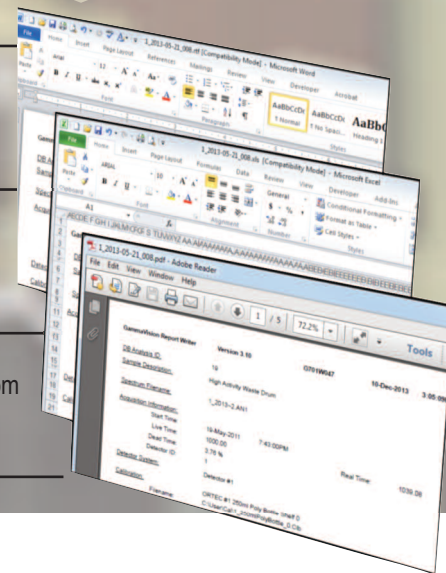

Version 3.10 Updates

- New!** 64-Bit Windows 7, Crystal Reports 11.5, and GammaVision 8 Compatibility.
- New!** Export compatibility with newer file formats (Word, Excel, Rich Text, and PDF).
- New!** Over 75 new Analysis Parameters Available.
- New!** Crystal Report 11.5 Templates with New Analysis Parameters.
- New!** Improved Printing and Report Export capability.

GammaVision Report Writer

Reporting Options

- γ Different Report Templates by Application or Sample Type
- γ Automatic Print Option
- γ Include Spectrum Image Report Option
- γ Spectrum Image Report Format:
 - Customizable Plot Title and Color Scheme
 - Optional Axis Labels and Sample Description
 - Acquisition Date/Time and Real/Live Time Reported
 - Regions of Interest Outlined or Filled
 - Horizontal Axis
 - Energy or Channels
 - Tick Marks or Grid Lines
 - Full Scale or User-Defined Range
 - Vertical Axis
 - Log or Linear Scale
 - Tick Marks or Grid Lines
 - Auto Scale or User-Defined Maximum
- γ Custom Reports
 - Example Crystal Report Templates are provided to produce lengthy detailed reports or a short summary of key data.
 - Extensive report customization available through Crystal Reports 11 or 11.5 as well as Visual Studio 2010.
 - Microsoft Access Reporting functions can be enabled with the Report Writer database. (Requires Microsoft Access)

Report Export Options

- γ Automatic Export on Analysis Option
- γ Export File Format:
 - Microsoft Word and Excel
 - PDF
 - ASCII Text
 - Rich Text
 - Crystal Reports
 - Additional formats available for manual export from the Crystal Reports Built-In Viewer Control.

Data Management

- γ Query Samples by Spectrum Name, Date Range, Sample Description and Detector by exact match or wild card criteria.
- γ Move/Copy selected analysis data and (optionally) associated files to a specified database and file folder for Archive functionality.
- γ Delete selected analysis data and (optionally) associated files for system Purge and Clean-up functionality.
- γ Use Microsoft Access or custom applications to run queries, combine multiple sample results, generate exception reports, or integrate with Laboratory Information Management systems.

Specifications

Operating System: Windows 7 (32 and 64-bit) and Windows XP (32-bit) Operating Systems.

Report Template Editors: Crystal Reports XI and Visual Studio 2010 (not included with Report Writer).

Database Format: Microsoft Access.

Ordering Information

Model	Description
A44-BW	Report Writer for GammaVision.
A44-UW	Report Writer Update from A44-B32 or A44-BW models.
A44-GW	Documentation for GammaVision Report Writer

Specifications subject to change
092718

ORTEC[®]

www.ortec-online.com

Tel. (865) 482-4411 • Fax (865) 483-0396 • ortec.info@ametek.com
801 South Illinois Ave., Oak Ridge, TN 37830 U.S.A.
For International Office Locations, Visit Our Website

AMETEK[®]