


Introducing GammaVision Report Writer V3.10

Robust GammaVision Analysis. . .

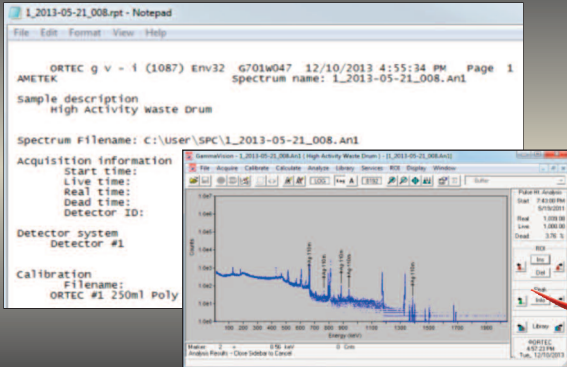

GammaVision Report Writer

Reporting Options

- γ Different Report Templates by Application or Sample Type
- γ Automatic Print Option
- γ Include Spectrum Image Report Option
- γ Spectrum Image Report Format:
 - Customizable Plot Title and Color Scheme
 - Optional Axis Labels and Sample Description
 - Acquisition Date/Time and Real/Live Time Reported
 - Regions of Interest Outlined or Filled
 - Horizontal Axis
 - Energy or Channels
 - Tick Marks or Grid Lines
 - Full Scale or User-Defined Range
 - Vertical Axis
 - Log or Linear Scale
 - Tick Marks or Grid Lines
 - Auto Scale or User-Defined Maximum
- γ Custom Reports
 - Example Crystal Report Templates are provided to produce lengthy detailed reports or a short summary of key data.
 - Extensive report customization available through Crystal Reports 11 or 11.5 as well as Visual Studio 2010.
 - Microsoft Access Reporting functions can be enabled with the Report Writer database. (Requires Microsoft Access)

Report Export Options

- γ Automatic Export on Analysis Option
- γ Export File Format:
 - Microsoft Word and Excel
 - PDF
 - ASCII Text
 - Rich Text
 - Crystal Reports
 - Additional formats available for manual export from the Crystal Reports Built-In Viewer Control.

Data Management

- γ Query Samples by Spectrum Name, Date Range, Sample Description and Detector by exact match or wild card criteria.
- γ Move/Copy selected analysis data and (optionally) associated files to a specified database and file folder for Archive functionality.
- γ Delete selected analysis data and (optionally) associated files for system Purge and Clean-up functionality.
- γ Use Microsoft Access or custom applications to run queries, combine multiple sample results, generate exception reports, or integrate with Laboratory Information Management systems.

Specifications

Operating System: Windows 7 (32 and 64-bit) and Windows XP (32-bit) Operating Systems.

Report Template Editors: Crystal Reports XI and Visual Studio 2010 (not included with Report Writer).

Database Format: Microsoft Access.

Ordering Information

Model	Description
A44-BW	Report Writer for GammaVision.
A44-UW	Report Writer Update from A44-B32 or A44-BW models.
A44-GW	Documentation for GammaVision Report Writer
Options	
A44-CR	Crystal Reports Professional

Specifications subject to change
121917