



Introducing GammaVision Report Writer V3.10

Robust GammaVision Analysis...

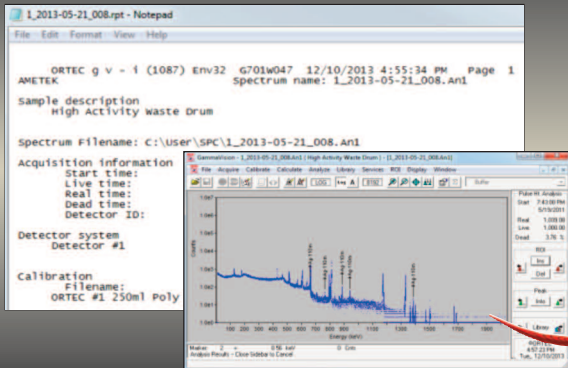

... in Brilliant Custom Reports

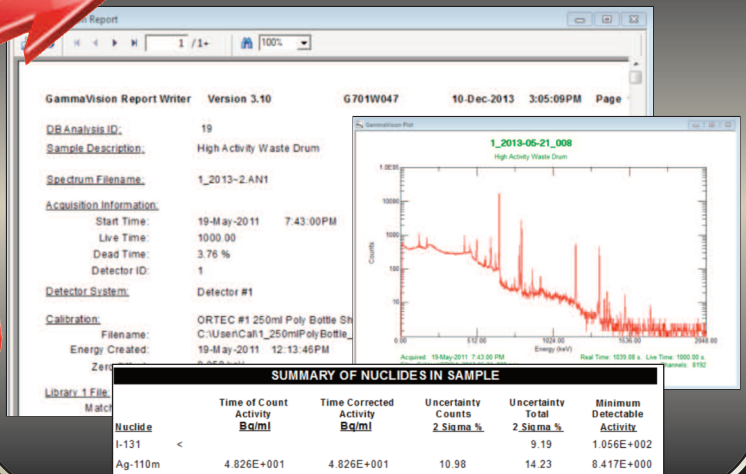

Export to Word, Excel, PDF, and More!

Why Report Writer?

Compatibility


- Windows 7 (32 & 64 bit) and Windows XP
- Microsoft Access Database
- Crystal Reports 11.5
- GammaVision 8 and earlier

Process Efficiency

- Automated Custom Reports with Spectrum Image
- Extensive Report Export Options (Word, Excel, PDF)
- Simplified LIMS Data Integration
- Data Management for Search, Archive, Purge functions

Defendable Results

- Populated directly from GammaVision Analysis
- Automated Calculations and Data Validation through Custom Report Templates
- Analysis Data supported by Spectrum Graphics

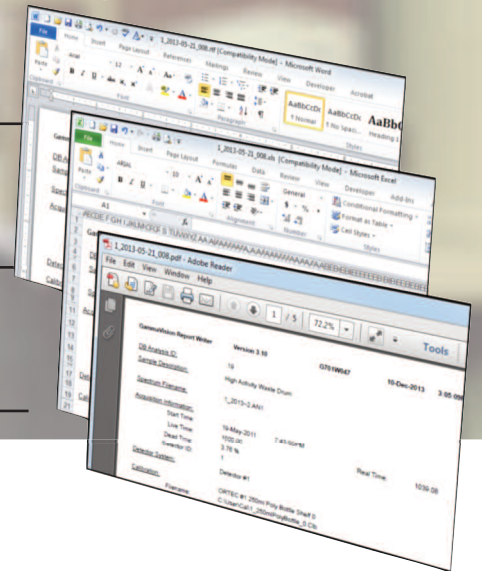

Version 3.10 Updates

New! 64-Bit Windows 7, Crystal Reports 11.5, and GammaVision 8 Compatibility.

New! Export compatibility with newer file formats (Word, Excel, Rich Text, and PDF).

New! Over 75 new Analysis Parameters Available.

New! Crystal Report 11.5 Templates with New Analysis Parameters.

New! Improved Printing and Report Export capability.

GammaVision Report Writer

Reporting Options

- γ Different Report Templates by Application or Sample Type
- γ Automatic Print Option
- γ Include Spectrum Image Report Option
- γ Spectrum Image Report Format:
 - Customizable Plot Title and Color Scheme
 - Optional Axis Labels and Sample Description
 - Acquisition Date/Time and Real/Live Time Reported
 - Regions of Interest Outlined or Filled
 - Horizontal Axis
 - Energy or Channels
 - Tick Marks or Grid Lines
 - Full Scale or User-Defined Range
 - Vertical Axis
 - Log or Linear Scale
 - Tick Marks or Grid Lines
 - Auto Scale or User-Defined Maximum
- γ Custom Reports
 - Example Crystal Report Templates are provided to produce lengthy detailed reports or a short summary of key data.
 - Extensive report customization available through Crystal Reports 11 or 11.5 as well as Visual Studio 2010.
 - Microsoft Access Reporting functions can be enabled with the Report Writer database. (Requires Microsoft Access)

Report Export Options

- γ Automatic Export on Analysis Option
- γ Export File Format:
 - Microsoft Word and Excel
 - PDF
 - ASCII Text
 - Rich Text
 - Crystal Reports
 - Additional formats available for manual export from the Crystal Reports Built-In Viewer Control.

Data Management

- γ Query Samples by Spectrum Name, Date Range, Sample Description and Detector by exact match or wild card criteria.
- γ Move/Copy selected analysis data and (optionally) associated files to a specified database and file folder for Archive functionality.
- γ Delete selected analysis data and (optionally) associated files for system Purge and Clean-up functionality.
- γ Use Microsoft Access or custom applications to run queries, combine multiple sample results, generate exception reports, or integrate with Laboratory Information Management systems.

Specifications

Operating System: Windows 7 (32 and 64-bit) and Windows XP (32-bit) Operating Systems.

Report Template Editors: Crystal Reports XI and Visual Studio 2010 (not included with Report Writer).

Database Format: Microsoft Access.

Ordering Information

| Model | Description |
|----------------|---|
| A44-BW | Report Writer for GammaVision. |
| A44-UW | Report Writer Update from A44-B32 or A44-BW models. |
| A44-GW | Documentation for GammaVision Report Writer |
| Options | |
| A44-CR | Crystal Reports Professional |

Specifications subject to change
121917

ORTEC[®]

www.ortec-online.com

Tel. (865) 482-4411 • Fax (865) 483-0396 • ortec.info@ametek.com
801 South Illinois Ave., Oak Ridge, TN 37830 U.S.A.
For International Office Locations, Visit Our Website

AMETEK[®]